



IPEX FROM NSK  
DIGITAL APEX LOCATOR

iPex

instant accurate measurements



The new NSK Advanced Digital Apex Locator has been developed to accurately measure the length of canals in all conditions, including dry, wet and where blood is present.



Powerful Partners®




## Precision Digital Apex Locator

NSK's new digital apex locator insures great accuracy in all canal conditions and eliminates problems apparent in other Apex Locators which use analogue signals. The NSK iPex multi-frequency measuring software technology filters the irrelevant signals that can appear in curved canals or under other extraordinary circumstances. iPex also automatically select the best possible combinations of frequency to meet patient's canal condition.

- Instant accurate measurement of canal length
- User friendly interface
- Large Hi-contrast LCD screen for easy monitoring
- Compact and smart design
- Audible warning system
- Auto power shut-off
- AAA battery operation
- No manual calibration necessary

## Large Display Plus Audio Signals for Immediate Monitoring

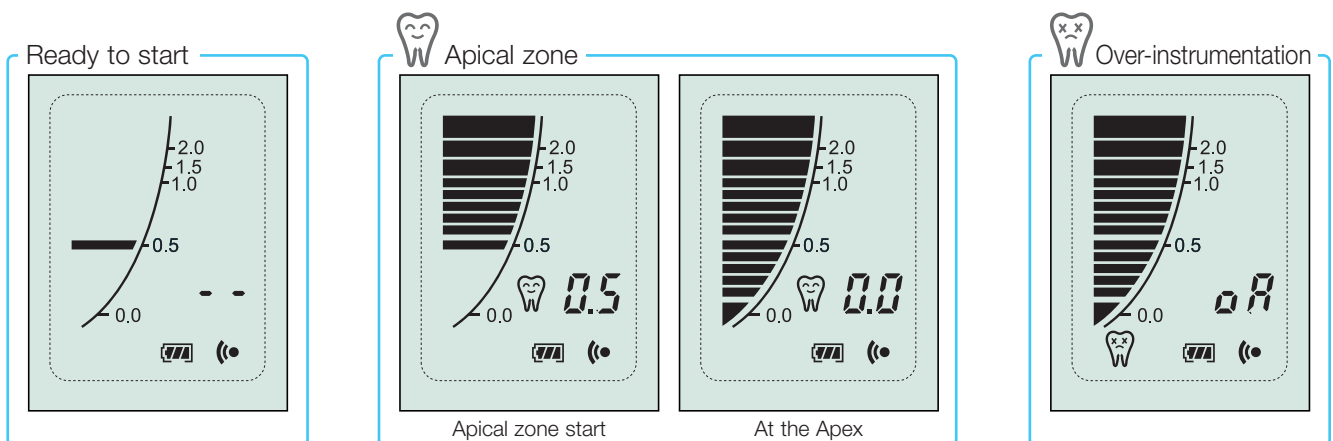

## iPex Complete Pack

MODEL ORDER CODE

iPex Y141-251

### Contents

- Control unit • Probe • File clip (3 pcs.) • Lip hook (3 pcs.)
- Extension cord • AAA manganese dry batteries (4 pcs.)

### Specifications

- Power Supply : for 1.5V AAA (LR03 or R03) batteries
- Continuous use : Approx. 50 hours • Dimensions : W80 x D80 x H95 (mm)
- Weight : Unit 185g



File clip, lip hook and extension cord can be autoclaved at 135°C max.

**NAKANISHI INC.**  
www.nsk-inc.com

700 Shimohinata Kanuma-shi  
Tochigi 322-8666,  
Japan  
TEL : +81 (0)289-64-3380  
FAX : +81 (0)289-62-5636

**NSK Europe GmbH** EC REP  
www.nsk-europe.de

Elly-Beinhorn-Strasse 8  
65760 Eschborn,  
Germany  
TEL : +49 (0)6196 77606-0  
FAX : +49 (0)6196 77606-29

**NSK United Kingdom Ltd**  
www.nsk-uk.com

Office 5, Gateway1000,  
Arlington Business Park, Whittle Way,  
Stevenage, SG1 2FP, UK  
TEL : +44 (0)1438 310670  
FAX : +44 (0)1438 310671

**NSK Oceania Pty Ltd**  
www.nsk-inc.com

Unit 22, 198-222 Young St.  
Waterloo, Sydney,  
NSW 2017, Australia  
TEL : +61 (0)2 8306 3000  
FAX : +61 (0)2 8306 3030

**NSK America-Canadian Branch**  
www.nsk-inc.com

1054 Centre St, Suite 283  
Thornhill, Ontario,  
Canada L4J 8E5  
TEL/FAX : +1 905 780 6757  
Toll Free : +1 866 780 6757

**NSK Middle East**  
www.nsk-inc.com

Room 6EA-701, 7th Floor, East Wing No.6  
Dubai Airport Free Zone,  
PO Box 54316 Dubai, UAE  
TEL : +971 (0)4-204-5507  
FAX : +971 (0)4-204-5508